

Cool Cube™ 08 at Frozen Temps (Varicella, MMRV, Zoster & FFP)

Model: FT-08

The Cool Cube™ 08 at Frozen Temperatures is not an ordinary transport cooler. It is an easy-to-carry, qualified container and pack-out system that keeps temperature sensitive product (varicella, MMRV, zoster, Fresh Frozen Plasma) frozen between -50° and -15°C for over 2.5 days – without using ice or electricity! This cooler is ideal for transport, but also excels at protecting valuable product during off-site use and power outages, regardless of hot conditions. The simple, validated, patented technology is the best compact thermal protection product available to maintain the cold chain!




4 Crêdo™ inside


- 1 Adjustable shoulder strap with pad is removable.
- 2 Double zipper closure for easy access.
- 3 Clear external pockets for documents and/or Digital Data Logger (DDL).
- 4 Six removable, temperature-controlling, PCM panels ensure a simple, quick & effective pack-out.

Qualified Pack-Out System


Remove PCM Panels

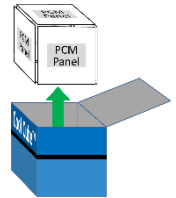


Prep PCM Panels

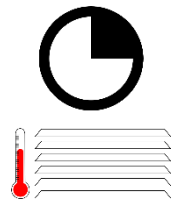

Assemble PCM Panels

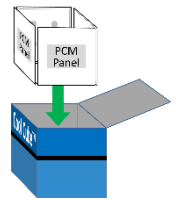

Pack-Out Product

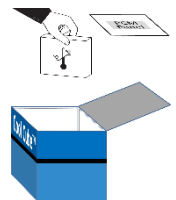

Close Bag



Features

- 6 patented PCM panels (Crêdo™ Series 20M) manage the temp
- Temp-Shield™ vacuum insulated panel (VIP) system
- Rugged, blue carrying case with a shoulder & hand strap
- Validated under ISTA Test Procedure 7D (Int'l Safe Transit Assoc)

Specifications

Temperature Range
-50° to -15°C (-58° to 5°F)

Performance
60+ Hours

Payload Area
11¼" x 6½" x 6½"

Exterior Dimensions & Weight
16" x 11½" x 11½" @ 20lbs

Go to Web Page >>



info@VeriCorMed.com
(866) 469-6019

For more information & pricing: www.VeriCorMed.com

