

# Cool Cube™ 08 at Room Temps (FFPE, platelets, biospecimens)

Model: RT-08



The Cool Cube™ 08 at Room Temperatures is not an ordinary transport cooler. It is an easy-to-carry, qualified container and pack-out system that keeps temperature sensitive product (FFPE, platelets, biospecimens) between 15° and 25°C for over 3 days – without using ice or electricity! This cooler is ideal for transport, but also excels at protecting valuable product during off-site use and power outages, regardless of hot or freezing conditions. The simple, validated, patented technology is the best compact thermal protection product available to maintain the cold chain!



4 Crēdo™ inside

## Features


- 6 patented PCM panels manage the temp
- Temp-Shield™ vacuum insulated panel (VIP) system
- Rugged, blue carrying case with a shoulder & hand strap
- Validated under ISTA Test Procedure 7D (Int'l Safe Transit Assoc)

## Specifications


- Temperature Range*  
15° to 25°C (59° to 77°F)
- Performance*  
83+ Hours (summer)
- Payload Area*  
11¼" x 6½" x 6½"
- Exterior Dimensions & Weight*  
16" x 11½" x 11½" @ 17lbs

## Validated Pack-Out System


Remove PCM Panels

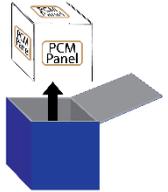


Prep PCM Panels

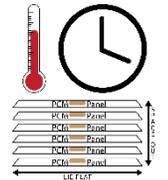

Assemble PCM Panels

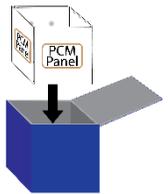

Pack-Out Product

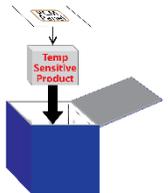

Close Bag



- 1 Adjustable shoulder strap with pad is removable.
- 2 Double zipper closure for easy access.
- 3 Clear external pockets for documents and/or Digital Data Logger (DDL).
- 4 Six removable, temperature-controlling, PCM panels ensure a simple, quick & effective pack-out.

Go to Web Page >>



info@VeriCorMed.com  
(866) 469-6019

For more information & pricing: [www.VeriCorMed.com](http://www.VeriCorMed.com)

