

Cool Cube™ 96 at Room Temps (FFPE, platelets, biospecimens)

Model: RT-96



The Cool Cube™ 96 at Room Temperatures is not an ordinary transport cooler. It is a lockable, wheeled, qualified container and pack-out system that keeps temperature sensitive product (FFPE, platelets, biospecimens) between 15 and 25 °C for over 5½ days (or between 20 and 24 °C for over 3½ days) – without using ice or electricity! This cooler is ideal for transport, but also excels at protecting valuable product during off-site use and power outages, regardless of hot or freezing conditions. The simple, validated, patented technology is the best compact thermal protection product available to maintain the cold chain!

Features

- 6 patented room PCM panels manage the temp
- Temp-Shield™ vacuum insulated panel (VIP) system
- Rugged, military spec case with large lifting handles & wheels
- Validated under ISTA Test Procedure 7D (Int'l Safe Transit Assoc.)

Specifications

Temperature Range / Duration
20 to 24 °C (68 to 75 °F) / 91+ hrs.
15 to 25 °C (59 to 77 °F) / 143+ hrs.
(summer profiles used for hrs.)

Payload Area
18" x 18" x 18"


Exterior Dimensions & Weight
26½" x 26½" x 28¼" @ 114 lbs.




- 1 Folding handle for ease of movement.
- 2 Holder and flexible grommet enables the use of an external temp monitor.
- 3 Multi-directional locking wheel system.
- 4 Six removable, temperature-controlling, PCM panels ensure a simple, quick & effective pack-out.

Validated Pack-Out System


Remove
PCM Panels

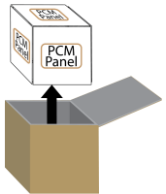


Prep
PCM Panels

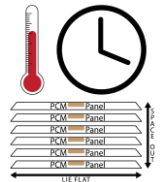

Assemble
PCM Panels

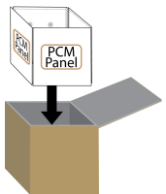

Pack-Out
Product

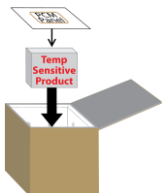

Close
Case



[Go to Web Page >>](#)



info@VeriCorMed.com
(608) 526-6901

For more information & pricing: www.VeriCorMed.com

